

Saint-Gobain Life Sciences -Bioprocess Solutions



Bio-Simplex® Media Bottle Assembly Systems

Closed-System Solutions for Aseptic Sample Handling, Fluid Transfer and Storage

Single-use technologies play a key role in reducing the risk of product cross contamination. Additionally, there are time and cost benefits associated with assembly, cleaning and validation. Bio-Simplex® Media Bottle Assembly Systems include a 2 port EZ Top® Container Closure, made from C-Flex®biopharmaceutical thermoplastic elastomer (TPE) resin and tubing, engineered to ensure an unobstructed fluid path with a secure elastomericseal between bottle and cap surfaces.

Designed for the bioprocess market, Media Bottle Assembly Systems are easily used with sterile connecting devices and thermal tubing sealers. Validation has been performed to demonstrate that the assemblies achieve a Sterility Assurance Level (SAL) of 10-6 per ANSI/AAMI/ISO 11137. Customer-generated designs, sizes, packaging and other options — such as connectors, heat-sealed ends or alternate materials — are available as custom orders.

Features and Benefits

- Closed system design that reduces contamination risk
- Low gas permeability
- · Low protein binding
- · Clarity for easy content viewing
- Volume graduated design
- Tubing welding / sealing capability sterile

Single-use Applications

- Media, sera, and buffers storage
- Media sterile addition (from Bio-Simplex Carboy to Media Bottle)
- Aseptic sampling




The Media Bottle Assembly System includes:

- Two port EZ Top® container closure
- C-Flex® tubing 18" (46cm), 1/8" (3.2mm) I.D., and 1/4" (6.4mm) O.D. with a plug end and dip tube for easy aseptic transfer
- 0.2µm vent filter for sterile gas exchange
- One media bottle (60mL to 2L)

BIO-SIMPLEX MEDIA BOTTLE ASSEMBLY SYSTEMS		
VOLUME	PART NUMBER	QTY / PACK
60mL	SA24MM-60-2F	10
125mL	SA38MM-125-2F-DB	10
250mL	SA38MM-250-2F-DB	10
500mL	SA38MM-500-2F-DB	10
1000mL (1L)	SA38MM-1L-2F-DB	5
2000mL (2L)	SA53B-2F-DB	5

Bio-Simplex Media Bottle Assembly Systems

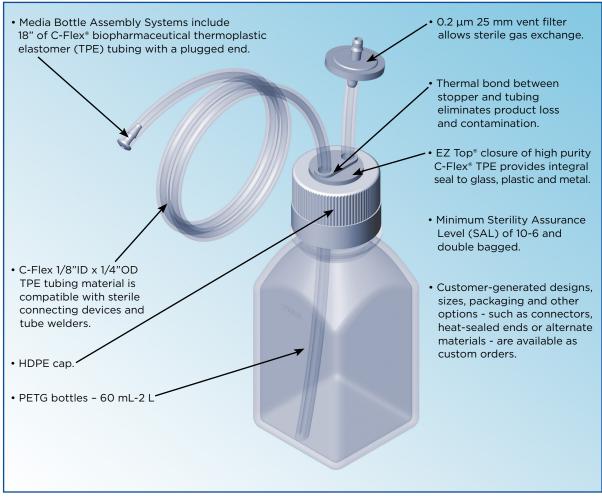

Bio-Simplex*, C-Flex* and EZ Top* are registered trademarks.



Saint-Gobain Performance Plastics

4451 110th Avenue North Clearwater, FL 33762 Tel: (727) 531-4191 Fax: (727 530-5603

www.biopharm.saint-gobain.com

<u>Contact us today for:</u> <u>Consultations • Samples • Quotes • Orders • On-Time Deliveries</u>

For registered access to Saint-Gobain Product Validation Summaries, CLICK HERE

IMPORTANT: It is the user's responsibility to ensure the suitability and safety of Saint-Gobain Performance Plastics products for all intended uses and that the materials to be used comply with all applicable medical regulatory requirements. Saint-Gobain Performance Plastics assumes no responsibility for any product failures that occur due to misuse of the materials it provides arising out of the design, fabrication or application of the products into which the materials are incorporated.

WARRANTY: For a period of 12 months from the date of first sale, Saint-Gobain Performance Plastics warrants this product to be free of defects in materials and workmanship. Our only obligation will be to replace any portion proving defective, or at our option, to refund the purchase price thereof.

SAINT-GOBAIN PERFORMANCE PLASTICS DISCLAIMS ALL IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE.